

# Club



**VERMEIREN**  
we care for you

# Club Family


**A series of beds combining over 65+ years of expertise in home care solutions.** The Club family from **VERMEIREN** is designed based on input from users, rental companies and health care professionals. A simple and safe four-section bed to provide the best comfort and care at home.

Adjust all Club beds electrically with their three Linak™ standard motors. Adjust the height according to the needs of either user or caregiver for care or transfers.

The (anti-) Trendelenburg position is standard included to go the extra mile for the user and bring additional health benefits when needed.



The controller is easy to use for both carer and user but can be locked if needed for additional safety.

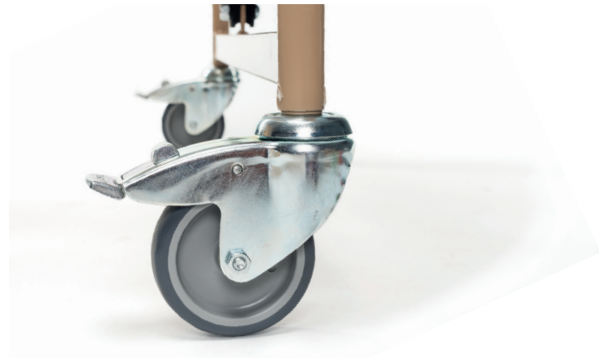


All 4 wheels can be locked for safety and security and allow for easy maneuvering.

## Club VARIO

**Club Vario is the height adjustable version of the Club bed.** On top of all the Club benefits, the base height can be set up in two different positions. Choose

the height that suits your situation or needs best, adjust electrically later to suit different needs throughout the day.

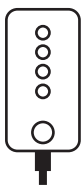

**Position 1: 33,5 cm** (electrically adjustable till 74,5 cm)

**Position 2: 40 cm** (electrically adjustable till 81 cm)

# Club

To ensure the best solution for you and your comfort, the Club bed has been designed with many configuration options. From the slatted base to the finishing covers.

**DISCOVER MORE** about Club's features below.



## Slatted base




Metal slats




Wooden slats

## Optional rastomat

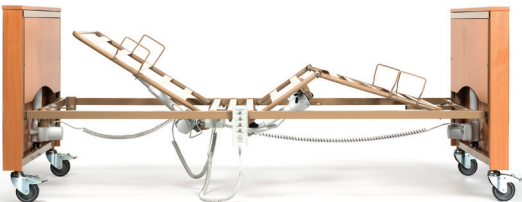

Without the rastomat

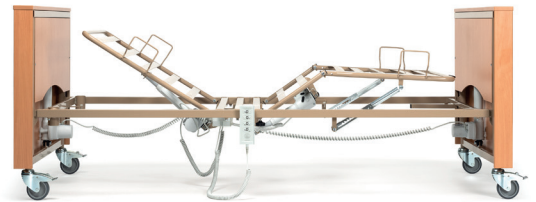

With the rastomat

## Side rails



Metal side rails



Split side rails (metal or wood)



Wooden side rails

## Wooden cover panels

Create a modern and striking look by adding wooden cover panels to Club. Perfect to create a cosier and clean look for your care bed.




## Options

We've selected the best **Club** options for you. Discover all options here.



A monkey pole adds additional support to the independency of patients when lifting or repositioning themselves.

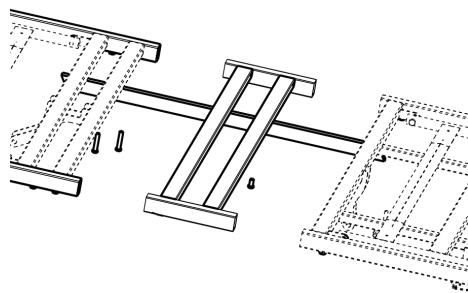

Club beds are created to suit everyone's needs! Easily add an additional length of 20 cm to the Club bed with our bed extension.



Serum holder



Transport kit

Configurations may differ depending on your region. Please consult our order forms for the most up-to-date information regarding all possibilities.

## Products designed to suit you

**VERMEIREN** always strives to come up with the best solution available to encompass your personal needs. As technology advances, so do we and as such we can offer you the most up to date products. This is the result of our continuous research and experience since 1957.

## Quality guaranteed

Regarding quality, **Club** complies with the highest norms, laid down by certified test institutes in most European countries.

## Technical data Club & Club Vario

■ Total width	1005; 1037 mm	■ Max. lower leg adjustment	-16,5°
■ Total height (without lifting pole)	834 - 1262 mm	■ Mattress spec.	RG 35, thickness 120 mm, 2000x900 mm
■ Total length	2125; 2175 mm	■ Max. load lifting pole	80 kg
■ Height adjustment of the bed surface (without mattress)*	400 - 810 mm	■ Weight	60,8 kg; 62 kg
■ Max. backrest adjustment	69°	■ Maximum user weight	165 kg
■ Max. thigh adjustment	20,7°	■ Additional info	CE

\* Club Vario pos.2. Measurements in mm and kg. Deviations in sizes and weight (15 mm and 1,5 kg) are possible. Subject to technical alterations. Images for illustration purposes only. Certain features may be optional.

## VERMEIREN Group

Vermeirenplein 1-15  
B-2920 Kalmthout  
Belgium

[www.vermeiren.com](http://www.vermeiren.com)  
[info@vermeiren.com](mailto:info@vermeiren.com)



For more information, contact your local dealer

**VERMEIREN**  
we care for you