



**“WHEELCHAIRS THAT MAKE
A GOOD FIT”**

“Roll towards your dreams!”

- VERMEIREN

100% PERSONALIZABLE WHEELCHAIRS

WHY TNS?

At TNS, we take pride in being the premier choice for those seeking top-quality, reliable, and knowledgeable products in the made-to-measure market. Our commitment to excellence is unmatched, and we strive to make a positive difference with every product we manufacture. When you choose TNS, you're not just selecting a handmade personalizable product – you're choosing a trusted partner that truly understands your unique needs and demands.

Each TNS wheelchair is made-to-measure, by hand, in the Netherlands, so we recommend an initial assessment with a local **VERMEIREN** sales representative. TNS works with fully trained specialists to ensure that all assessments are carried out before producing the wheelchair that fits your needs, and your wishes.




THE BEAUTY OF PERSONALIZATION

Experience the ultimate in personalization with our wheelchair fitting service. Not only will your chair fit you perfectly, but you can also choose from a wide range of colours and customizable details to make your chair uniquely your own. Our goal is to create a chair that truly represents your style and personality, making every moment of your day as comfortable and stylish as possible.



Lightweight

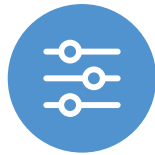

100%
Personalizable



Handmade



HANDMADE IN THE **NETHERLANDS**

"TNS is a company specialized in making amazing wheelchairs and listens to their customers wishes from the start. Nothing is Impossible! You do not feel that you are just a "number" and they are searching and listening for the best options for you and your new wheelchair. I needed a wheelchair that is small and easy to handle but also strong because I travel a lot for my sport (table tennis). Till today, I am still grateful they are supporting me on my way to the Paris 2024 Paralympic Games."

-Kelly





DISCOVER THE COLLECTION



Unlimited possibilities

**Activator
P7**

**Activator 2.5
P8**

**Activator F.C.
P9**

**Proval
P10**

**RSX
P11**

Activator

The Open Frame



Open frame

30 mm tube

5 kilograms

transport weight

Minimal-frame / L-frame

Adjustable seat mat

Tension-adjustable backrest with breathable fabric

Push/pull or scissor brakes

Choice of 200 frame colors in matte, gloss or glitter

Open-frame active wheelchairs are usually lighter in weight compared to closed-frame wheelchairs, making them easier to maneuver and transport. Activator allows for a lower center of gravity, which provides greater stability and control when navigating rough or uneven terrain.



Activator 2.5

The Most Rigid



Box-frame wheelchairs have a rigid, welded frame that forms a box-like structure around your body. The Activator 2.5 design provides a direct transfer of energy from the user's upper body to the wheels, which can result in better performance and efficiency. This can be particularly important for users who rely heavily on their upper body strength to propel themselves.



Box frame

25 mm tube

4,5 kilograms

transport weight



Activator F.C.

The Most Foldable

Foldable front frame

30 mm tube

5,5 kilograms

transport weight



The ability to fold the wheelchair allows for greater convenience and flexibility in terms of storage and transportation. Activator F.C. is the perfect solution if you require a high degree of portability and convenience, but still want the benefits of an active wheelchair.

Proval

The Lightest



Box frame

15 x 30 mm tube

4 kilograms

transport weight



Proval design combines all benefits from a box frame, with a unique oval frame tube shape that makes it the lightest wheelchair in our range. This unique feature also provides Proval with a minimalist and discreet design where your lifestyle comes first.



RSX

The Most Sporty look

Rear Suspension Experience

25 mm tube

6,5 kilograms

transport weight



Our RSX (Rear Suspension Experience) active wheelchair design is the epitome of sportiness. With its front shape frame, you can enjoy additional freedom without any dominant frame tubes limiting your movement or impeding your transfers. Additionally, the RSX offers a rear suspension experience that enhances performance and ensures a more comfortable ride.



TABLE OF COMPARISON

Frame

Transport weight from
(without wheels and cushion)

Foldable backrest

Construction

Crash tested

ACTIVATOR

The Open Frame



30 mm tube

Open frame

5 kg



Aluminium



ACTIVATOR 2.5

The Most Rigid



25 mm tube







Box frame


4,5 kg



Aluminium



ACTIVATOR F.C.	PROVAL	RSX
The Most Foldable	The Lightest	The Most Sporty
 <p>30 mm tube</p> <p>Foldable front frame</p> <p>5,5 kg</p> <p>✓</p> <p>Aluminium</p> 	 <p>15 x 30 mm tube</p> <p>Box frame</p> <p>4 kg</p> <p>✓</p> <p>Aluminium</p> 	 <p>25 mm tube</p> <p>Suspension</p> <p>6,5 kg</p> <p>✓</p> <p>Aluminium</p> 

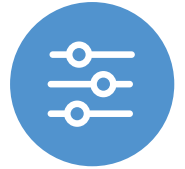



Activator 2.5

DYNAMIC

70000k Role



Baton

CRITERIA WHEN CHOOSING AN ACTIVE WHEELCHAIR



Activity

The first thing to consider is your activity level. An active wheelchair is typically designed for individuals who lead an active lifestyle and require a high degree of mobility and functionality on difficult parcours.



Mobility needs

Your specific mobility needs should also be taken into account. For example, if you need to navigate rough terrain or climb stairs, a wheelchair with larger wheels or additional support features may be necessary.



Body measurements

The wheelchair should be sized to fit your body, including height, weight, and the shape of your body. This ensures maximum comfort and support while using the chair.



Lifestyle & preferences

Your lifestyle and personal preferences should also be considered. This includes factors such as preferred color, design, and features, as well as any specific accessories or personalized options that may be required.



**WE BUILD 100% PERSONALIZABLE
SPORTS WHEELCHAIRS FOR THE
BEST SPORT PERFORMANCE
YOU CAN GET!**





SPORTS WHEELCHAIRS



Our team has been crafting and customizing wheelchairs for several decades, and we're dedicated to meeting your specific needs for sports and other activities. By first listening to your requirements, we can create a wheelchair that enhances your performance in dance, tennis, basketball, and more. With a TNS wheelchair, you can focus on your activity without any worries, as our product becomes a valuable asset to help you excel.



Contact your sales representative to discover all possibilities.

VERMEIREN
we care for you



distributed by
VERMEIREN

NV VERMEIREN
Vermeirenplein 1-15 ,
2920 Kalmthout - Belgium
info@vermeiren.com